Upute za upotrebu

BEEABeauty

TB-2187

ALL-IN-ONE
uređaj za njegu lica
CANC, RF, LED, ION, Hot&Cold
ANTI AGING

UPOZORENJA

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Datum kupnje: _____ Pečat i potpis: _____

WARNING

-  This device is not a toy and is not suitable for children. Please keep it out of childrens reach.
-  Always verify the device before you start using it, do not use it if the device is damaged or broken.
-  If used dry (without face product) this device may generate static electricity. This is normal and does not affect its use.
-  Keep the USB charging cable and the device's charging station dry.

PRECAUTIONS

-  Do not use this device in the following situations:
 -  Pregnant women or while nursing; children; people who have a pacemaker or other electronic medical equipment; people with allergic reactions.
 -  If you have skin conditions, such as open wounds, ringworm, eczema, rash, burns, sunburns, abnormal wound healing, with makeup on etc.

- People with diabetes, lupus erythematosus, porphyria or congestive heart disease, bleeding disorder, history of immunosuppressive diseases (including HIV infection and AIDS), etc.
 - Do not use this device in conjunction with any medical electronic equipment such as a pacemaker, internal defibrillator, etc. Do not use it if you have an arrhythmia or any known heart disease.
 - Do not use this device if you have metal implants, piercings, jewelry, or moles in the area of use.
 - Do not use this device if you are receiving radiation therapy or chemotherapy treatments.
 - People with heating-stimulating diseases, such as herpes simplex virus.
- ⚠ Do not use this device if you take any of the following medications:
- If your skin is currently being treated with or has recently been treated with Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), topical isotretinoin, and azelaic acid.
 - If you have taken any form of Accutane A or Roaccutane A within the past six months.
 - If you are taking photosensitizing agents or any other doctor-prescribed medication, please check the patient information leaflet before you start using this device. If the patient information leaflet states that it can cause photo-allergic or photo-toxic reactions, or if it specifies avoiding sun exposure when taking the medicine, please do not use the device.

- If you take anticoagulant medication, including large doses of aspirin, only use the device one week after you stopped the treatment. Use the device exclusively according to the recommendations, avoiding botox-treated areas within 2-3 weeks after the botox treatment, and consult a dermatologist before use.

This device is suitable to use on the face (forehead, cheeks, chin). Never use it on the eyelids, lips, chest, genital area, throat, or on areas where you have had surgery or epilation. If your skin is light-sensitive or prone to allergic reactions, do not use light therapy. Do not stare at the light therapy source. Do not use light therapy on your skin if the treatment area has infections, eczema, burns, inflammation of hair follicles, open lacerations, herpes, wounds or lesions, hematomas, etc. Do not look directly at the LED light. Use eye protection during the light therapy treatment, such as sunglasses or protective goggles (not included). Do not repair, disassemble or modify the device by yourself, as improper handling may damage it. Do not use this device for any other purpose. Clean this device surface and accessories regularly to extend their lifespan. It is recommended to use a gel compatible with this device (purchase separately). Using any other cream or gel may affect the normal operation of this device and its RF function.

Device Description

TB-2187

🔍 Device Specifications

Name: All-In-One Facial care device BEEA Beauty, CANC, RF - Anti-Aging

Model No.: TB-2187

Input: 5V / 2A

Power: 3.7V / 8.5W

Type of Battery: Li-poly 3.7V 1500mAh

Device Components

① Main Unit

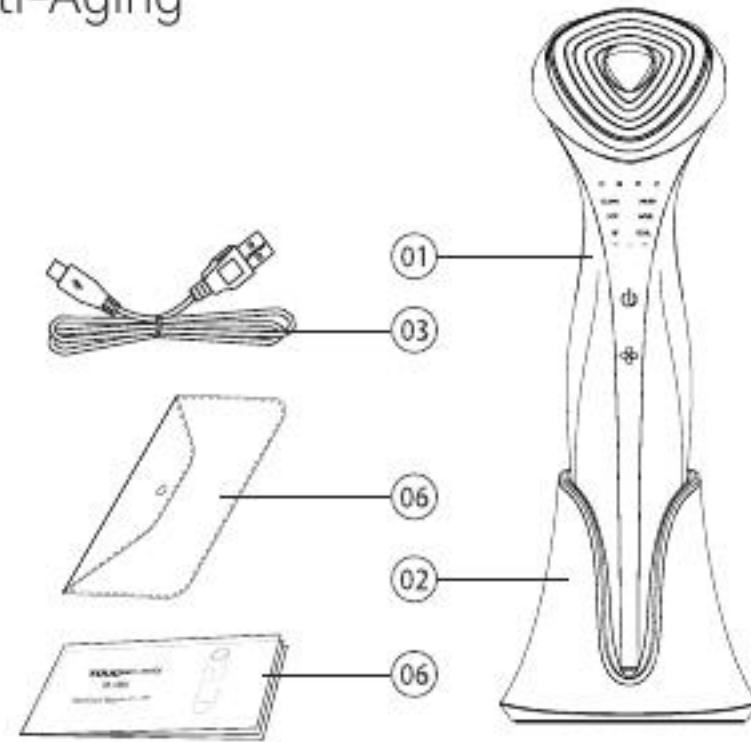
② Charging Dock

③ Charging Cable

④ Cleaning Cloth

⑤ Storage Bag


⑥ Instruction Manual

Picture 1

🔦 Device Features

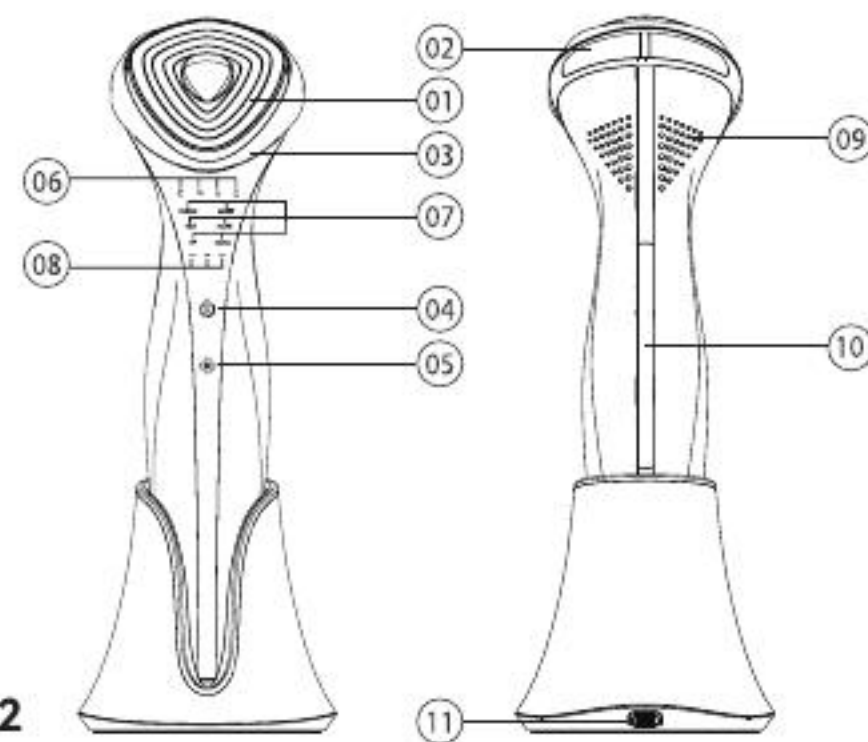- ① Facial Massage Head
- ② Partial Massage Head
- ③ Magnetic Cotton Ring
- ④ Power/Mode Switch Button
- ⑤ Intensity Level Button
- ⑥ CANC Indicator
- ⑦ Mode Indicator
- ⑧ Intensity Level Indicators
- ⑨ Cooling Vents
- ⑩ Electrode Sensor
- ⑪ Type-C Charging Terminal

Picture 2

Preparing for Use

TB-2187

※ The battery is low before you use the device for the first time. Charge the device with the USB cable.

🔋 Charging the Device

- This device features a wireless charging system. Before you start charging the device, insert one end of the charging cable into the charging dock and plug the other end into an adapter (not included). Connect the adapter to a wall socket. Ensure that you use an adapter that has an output current superior to 5V/1A.

Place the device on the charging base to start charging. (Picture 3)

When charging, the intensity level indicators flash in a cycle. When fully charged, all the intensity level indicators light up continuously. (Picture 4). Charging the device takes about 5 hours. A full battery allows for approximately 60 minutes of continuous use.


Picture 3

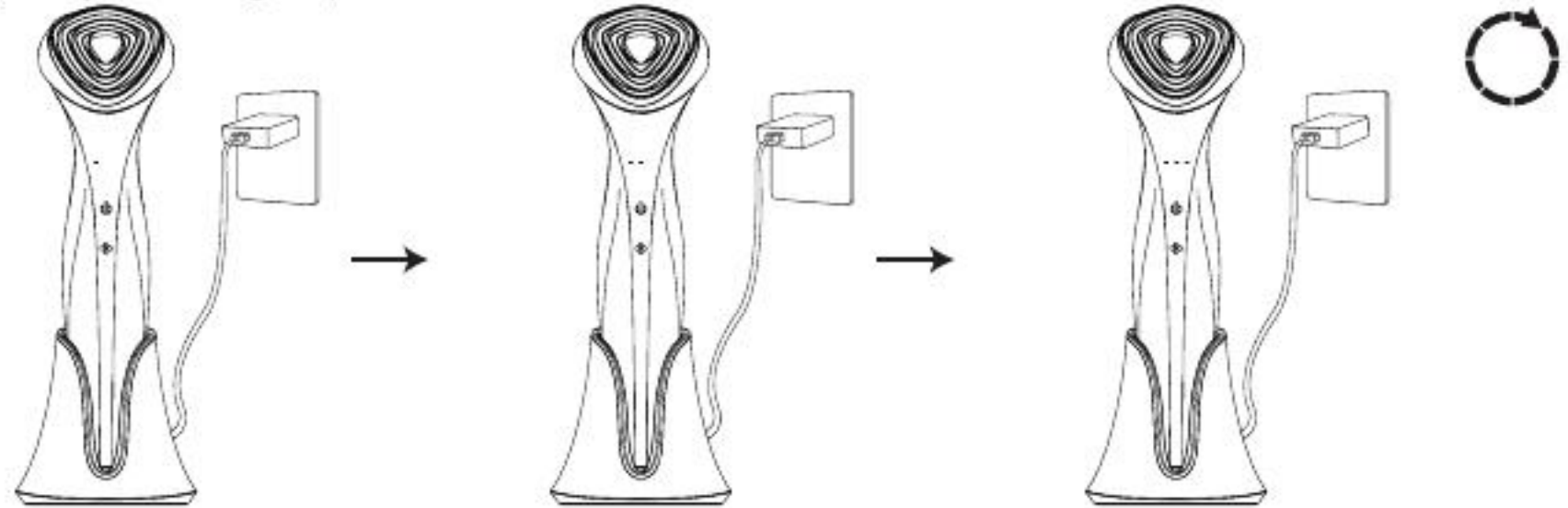- If the battery is low, the mode and intensity level indicators flash 5 times at the same time, and the device automatically shuts off.

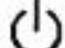
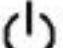
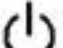
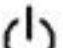
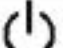
🔌 Turn ON-OFF & Mode Switch

- ✂ This device's power/mode switch  and intensity level buttons  have a touch-pressing design.
- When you press the button  with your finger, the device will vibrate to indicate that your finger is pressing the right area. You will need to apply pressure to turn it ON. If you apply pressure on the button and there is no vibration, you will need to re-adjust your finger position.
- The button  does not vibrate when you press it. Apply enough pressure to use it.


Picture 4

※ Operation Instructions


| Function | Function Introduction | Using Method | Intensity Level |
|-------------|--|--|---|
| Turn ON-OFF | | Press & hold  | \ |
| CLEAN Mode | Positive ions + vibrations | Turn ON the device for the first time CLEAN - default mode | Press  to adjust the intensity level - 3 levels. Press twice in quick succession to turn ON-OFF the vibration. |
| HOT Mode | EMS + hot therapy + vibrations | Press  1 time | Press  to adjust the EMS intensity level - 3 levels. Press twice in quick succession to turn ON-OFF the vibration. |
| RF Mode | RF + EMS + MC + red + infrared light | Press  1 time | Press  to adjust the RF intensity level - 3 levels. |
| MOIST Mode | Negative ions + vibrations | Press  1 time | Quickly press  to adjust the vibration intensity level - 3 levels. Press twice in quick succession to turn ON-OFF the vibration. |
| ACNE Mode | Red & blue mixed light + infrared light + hot massage | Press  1 time | \ |
| COOL Mode | Cold massage | Press  1 time | \ |

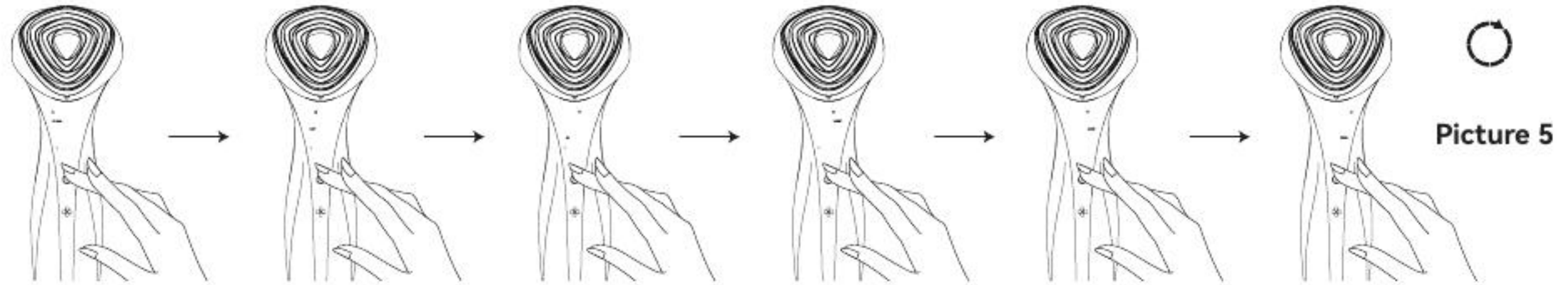- You can switch modes in a cycle: **CLEAN→HOT→RF→MOIST→ACNE→COOL**. When switching between modes, the light indicator of each corresponding mode lights up. (Picture 5)

- The intensity level can be adjusted in a cycle: **Level 1 → Level 2 → Level 3**. When you switch between intensity levels, the corresponding indicator lights up. If the intensity level of the mode you chose cannot be adjusted, the intensity level indicator turns OFF. (Picture 6)

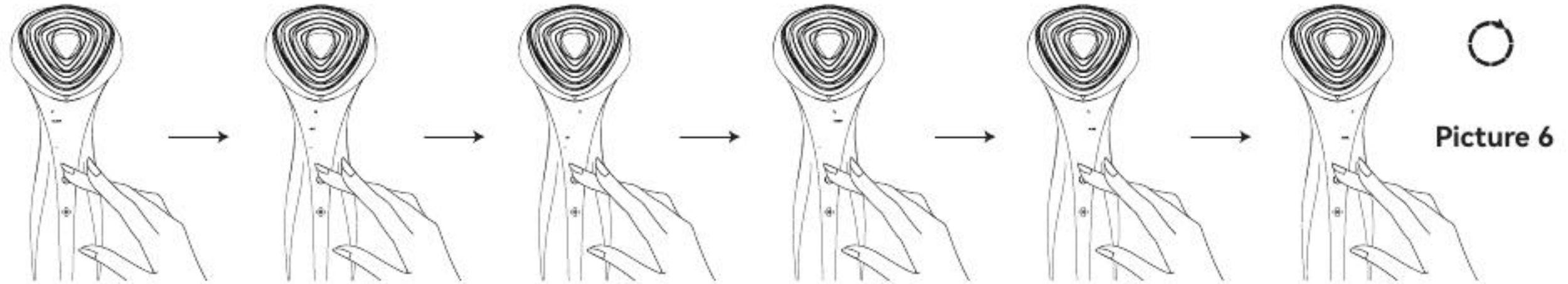👁️ Some functions of the partial massage head activate at the same time as the facial massage head. You can use the partial massage head according to your needs:

| Mode | CLEAN | HOT | RF | MOIST | ACNE | COOL |
|--------------------|----------------------------|----------------|-----------|---------------------------|------|------|
| Function Activated | Positive ions + vibrations | EMS+vibrations | RF+EMS+MC | Negative ions +vibrations | / | / |


⚠ ATTENTION

- When you use the partial massage head, ensure that the middle area is in contact with the skin so that both electrodes are in contact with the skin. (Picture 7)
- When you need to use the RF mode of the topical care head, use the lowest intensity level because the skin around the eyes is sensitive. Do not touch the eyelids and eyeballs when using this device.

※ When you use CLEAN/MOIST modes, ensure that your fingers are touching the electrode sensor while using the device. (Picture 8)

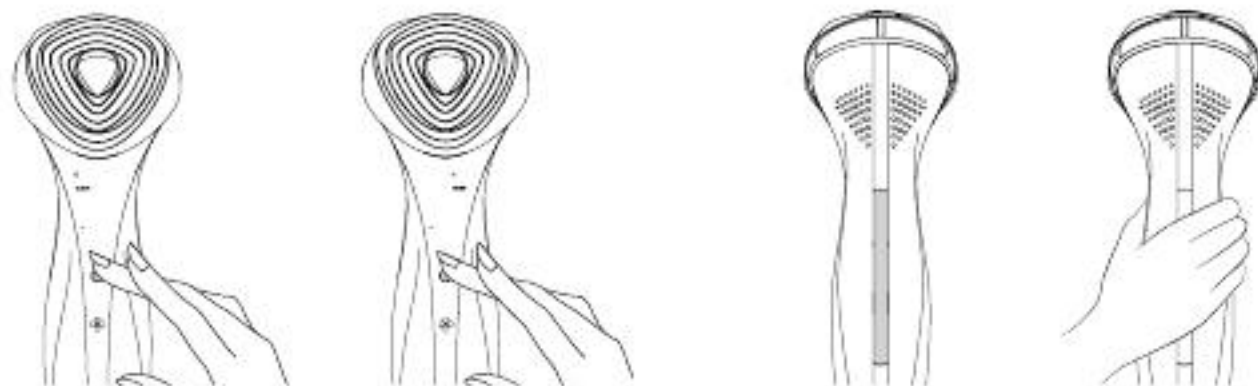Picture 7

Picture 8

※ This device has a memory function. It automatically activates the mode you have used before turning OFF the device. However, if the device runs out of battery, it activates the default mode setting after charging it.

🌙 Introducing CANC Scientific Skincare:

※ This mode combines the use of skincare devices and topical skincare products to boost efficiency. This 4-step scientific method Clean-Activate-Nourish-Care gives life to CANC skincare model through this device's functions:

※ This device does not include skincare products. You can use your favorite skincare products.

| CANC Clean – Activate – Nourish – Care 4-Step Skincare | Function & Recommended Frequency | Application with Skincare Products |
|--|-------------------------------------|--|
| 1.[C] Clean Deeply cleanses the skin. | CLEAN Mode (daily use) | Remove makeup, double-cleansing. Use with cotton pads and makeup remover. |
| 2.[A] Activate Activates the skin and opens pores. | HOT Mode (daily use) | Open pores and boost absorption. Use with toner or sheet mask. |
| 3.[N] Nourish Photoelectric beauty promotes absorption. | RF Mode (2-3 times per week) | Anti-aging: helps to firm the skin, use with a special gel. |
| | MOIST Mode (daily use) | Boosts absorption. After applying essence. |
| | ACNE Mode (as needed) | Reduces acne breakouts. Use directly on the skin or with acne skincare products. |
| 4.[C] Care Lock in nutrients and boosts hydration levels. | COOL Mode (daily use) | Helps to tighten pores, calms the skin. Use after applying lotion or moisturizer. |

How to Use

TB-2187

※ To ensure the correct operation of this device, please read the instructions carefully before use.

👁️ Using the Device

● Cleanse your skin before you use the device, and ensure that the device is charged.

※ Using Method:

| CANC | Mode | Function | Recommended Using Method | Using Technique |
|----------|-------|---|---|--|
| C | CLEAN | Cleansing Mode Positive ions +vibrations cleansing | <ol style="list-style-type: none">① Remove the magnetic ring. Place the cotton pad on the massage head surface and press the magnetic ring to fix the cotton pad.② Apply your regular makeup remover onto the cotton pad.③ Turn ON the device and choose CLEAN mode.④ Adjust the vibration intensity level according to your preference.⑤ Cleanse your skin with circular motions as illustrated.⑥ After use, turn OFF the device. |  |

| CANC | Mode | Function | Recommended Using Method | Using Technique |
|------|------|---|---|--|
| Q | HOT | Revitalize Mode: EMS + hot massage + vibrations | <ol style="list-style-type: none"> ① Apply toner or a sheet mask to your face. ② Turn ON the device and activate the HOT mode. ③ Adjust the EMS intensity level according to your skin tolerance. ④ Apply the facial massage head to your skin and massage as illustrated. Use this mode for 4-6 minutes/session. ⑤ After use, turn OFF the device. |  |
| U | RF | Anti-Aging Mode: RF(2MHz) + EMS + MC+ red light + infrared light | <ol style="list-style-type: none"> ① Evenly apply water-soluble gel to your face. ② Turn ON the device and activate the RF mode. ③ Adjust the RF intensity level according to your skin tolerance. ④ Massage the skin with the facial massage head in circular motions as illustrated. Use this mode for 4-6 minutes/session. ⑤ After use, wash your face. |  |

| CANC | Mode | Function | Recommended Using Method | Using Technique |
|------|-------|---|---|--|
| | MOIST | Boosting Mode: Negative ions + vibrations | <ol style="list-style-type: none"> ① Apply essence to the skin. ② Turn ON the device and activate the MOIST mode. ③ Adjust the vibration intensity level. ④ Massage your face with circular motions as illustrated. Use this mode for 4-6 minutes/session. ⑤ After use, turn OFF the device. |  |
| | ACNE | ACNE Mode: Red & blue mixed light + infrared light + hot massage | <ol style="list-style-type: none"> ① You can apply acne topical skincare products to acne-prone areas (not necessary, depends on your skin needs). ② Turn ON the device and activate ACNE mode. ③ Massage your face with slow circular motions as illustrated. Use it twice a day, for 5-min/session. ④ After use, turn OFF the device. |  |

| CANC | Mode | Function | Recommended Using Method | Using Technique |
|----------|------|-----------------------------------|---|---|
| C | COOL | COOL Mode: Cool therapy | ① Apply lotion or moisturizing cream to your skin. ② Turn ON the device and activate the COOL mode. ③ Apply the facial massage head to the skin and massage as illustrated. Use this mode for 4-6 minutes/session. ④ After use, turn OFF the device. |  |

Combining Modes:

- Daily Skincare: ① Wash face → ② CLEAN → ③ HOT → ④ MOIST → ⑤ COOL.
- Anti-Aging: It is recommended to use the RF 2-3 times a week. When using it, follow this recommended order ① wash face → ② CLEAN → ③ HOT → ④ RF → ⑤ MOIST → ⑥ COOL.

Use this device consistently for optimal results.

- Acne Breakouts: Do not use this mode if you do not have acne. If you have acne-prone skin or active breakouts, use this mode twice a day. Use ACNE mode before COOL mode.

⚠ ATTENTION

※ When using the radio-frequency function (RF mode):

- Remove all jewelry from your body before using the device, including earrings, rings, necklaces, etc.
- Re-apply gel if your skin gets dry during use.
- For first-time users, test the gel or essence on the skin 24h prior to initial use. Apply a small amount of gel or essence on the inside of the arm to ensure it will not cause a skin reaction.
- Use the device with the gel or the essence on the inside of the arm 24h before you use the device. After use, the skin may show signs of redness or appear slightly swollen for a short time. If these reactions do not disappear within 24 hours, please consult a doctor.
- Avoid using this device for too long on the cheeks or on areas you do not want volume reduction.
- RF function provides a warm sensation on the skin. If you feel the temperature is not enough for your skin, slow down the massage movements.
- During use, do not leave the device in one place without movement for longer than 2-3 seconds.
- Be sure to move the device across your skin to prevent minor burns due to different skin sensitivity and tolerance.
- Skin sensation and results may vary from person to person and area of use.

※ When using the light therapy function (RF/ACNE mode):

- Do not stare or look directly at the light source.

Cleaning & Maintenance

- ✕ Clean this device with a clean, slightly damp cloth before and after using the device.
- Always make sure that you turn OFF the device before cleaning.
- Do not immerse or clean the device directly with water or other liquid substances.
- Do not use steel wool, abrasive cleaning products, chemicals, or corrosive detergents to clean the device.

Storage

- This device has a built-in lithium battery. Do not place or store the device in hot and humid places. If you don't use the device frequently, make sure to charge the battery before you store it.

Disposal Information

- Lithium batteries contain substances that may pollute the environment. Do not dispose of this device with domestic waste. Instead, hand it over to an official recycling point to contribute to environmental protection.

 **Troubleshooting**

| Problem | Possible Cause | Solution |
|--|--|---|
| The device is not working. | The battery is running low or out of power. | Charge the device as guided in this instruction manual. |
| Cannot feel heat or microcurrent on the skin. | Not enough gel or essence was applied to the skin. | Apply more gel or essence. |
| | You are moving the treatment head too quickly. | Move the massage head slowly. |
| | The device stopped working. | Move the massage head slowly and evenly. |
| Note: If the device still does not work properly after using the above solutions, please stop using the device and contact the store or local distributor. | | |